# **Rat IL-4 ELISA Kit**





For the quantitative determination of rat interleukin 4(IL-4) concentrations in cell culture supernates, serum, and plasma.

This instruction must be read in its entirety before using this product.

For research use only. Not for use in diagnostic procedures.

# **Contact information**

Tel: 1-301-446-2499 Fax: 1-301-446-2413

Email: <u>tech@signalwayantibody.com</u> Web: <u>www.sabbiotech.com</u>

#### PRINCIPLE OF THE ASSAY

This assay employs the quantitative sandwich enzyme immunoassay technique. A monoclonal antibody specific for IL-4 has been pre-coated onto a microplate. Standards and samples are pipetted into the wells and any IL-4 present is bound by the immobilized antibody. Following incubation unbound samples are removed during a wash step, and then a detection antibody specific for IL-4 is added to the wells and binds to the combination of capture antibody-IL-4 in sample. Following a wash to remove any unbound combination, and enzyme conjugate is added to the wells. Following incubation and wash steps a substrate is added. A coloured product is formed in proportion to the amount of IL-4 present in the sample. The reaction is terminated by addition of acid and absorbance is measured at 450nm. A standard curve is prepared from seven IL-4 standard dilutions and IL-4 sample concentration determined.

### **DETECTION RANGE**

15.6 pg/ml - 1000 pg/ml

# **SENSITIVITY**

The minimum detectable dose was 7pg/mL.

# **SPECIFICITY**

This assay recognizes both natural and recombinant rat IL-4. The factors listed below were prepared at 50ng/ml in Standard /sample Diluent and assayed for cross-reactivity and no significant cross-reactivity or interference was observed.

Table 1: Factors assayed for cross-reactivity

Recombinant human	Recombinant mouse	Recombinant porcine	
IL-4	IL-4	IL-1α	
IL-4sR		IL-1β	
		IL-2	
		IL-6	
		TNF-α	

# **PRECISION**

The coefficient of variation of both intra-assay and inter-assay were less than 10%.

#### MATERIALS PROVIDED

- Aluminium pouches with a Microwell Plate coated with monoclonal antibody to rat IL-4 (8x12)
- 2 vials rat IL-4 Standard lyophilized,1000 pg/ml upon reconstitution
- 3. 2 vials concentrated Biotin-Conjugate anti-rat IL-4 antibody
- 4. 2 vials Streptavidin-HRP solution,
- 5. 1 bottle Standard /sample Diluent
- 6. 1 bottle Biotin-Conjugate antibody Diluent
- 7. 1 bottle Streptavidin-HRP Diluent
- 8. 1 bottle Wash Buffer Concentrate 20x (PBS with 1% Tween-20)
- 1 vial Substrate Solution
- 10. 1 vial Stop Solution
- 11. 3 pieces Adhesive Films
- 12. package insert

# **STORAGE**

Table 2: Storage of the kit

Unananad Kit	Store at 2 - 8° C. Do not use past kit expiration				
Unopened Kit	date.				
Opened/	Standard /sample Diluent				
	Concentrated Biotin-				
	Conjugate				
	Streptavidin-HRP solution	May be stored for up			
	Biotin-Conjugate antibody				
	Diluent	to 1 month at 2 - 8° C.**			
	Streptavidin-HRP Diluent				
	Wash Buffer Concentrate				
	20x				
	Substrate Solution				
	Stop Solution				
Reconstituted	Standard	Aliquot and store for up to			
Reagents		1 month at ≤20°C.			
		Avoid repeated freeze-			
	Otariaara	thaw cycles. Diluted			
		standard shall not be			
		reused.			
		Return unused wells to			
		the foil pouch containing			
		the desiccant pack,			
		reseal along entire edge			
		of zip-seal. May be stored			
		for up to 1 month at 2 - 8°			
		C.**			

\*\*Provided this is within the expiration date of the kit.

# THE REQUIRED ITEMS (not provided, but can help to buy):

- 1. Microplate reader (450nm).
- 2. Micro-pipette and tips: 0.5-10, 2-20, 20-200, 200-1000ul.
- 37 <sup>°</sup>C incubator, double-distilled water or deionized water, coordinate paper, graduated cylinder.

#### PRECAUTIONS FOR USE

- Store kit regents between 2<sup>°</sup>Cand 8<sup>°</sup>C. After use all reagents should be immediately returned to cold storage(2<sup>°</sup>C to 8<sup>°</sup>C).
- Please perform simple centrifugation to collect the liquid before use.
- 3. To avoid cross contamination, please use disposable pipette tips.
- 4. The Stop Solution suggested for use with this kit is an acid solution. Wear eye, hand, face, and clothing protection when using this material. Avoid contact of skin or mucous membranes with kit reagents or specimens. In the case of contact with skin or eyes wash immediately with water.

- Use clean, dedicated reagent trays for dispensing the washing liquid, conjugate and substrate reagent. Mix all reagents and samples well before use.
- 6. After washing microtiter plate should be fully pat dried. Do not use absorbent paper directly into the enzyme reaction wells.
- Do not mix or substitute reagents with those from other lots or other sources. Do not use kit reagents beyond expiration date on label.
- Each sample, standard, blank and optional control samples should be assayed in duplicate or triplicate.
  Adequate mixing is very important for good result. Use a minivortexer at the lowest frequency or Shake by hand at 10min interval when there is no vortexer.
- 9. Avoid microtiter plates drying during the operation.
- 10. Dilute samples at the appropriate multiple, and make the sample values fall within the standard curve. If samples generate values higher than the highest standard, dilute the samples and repeat the assay.
- 11. Any variation in standard diluent, operator, pipetting technique, washing technique, incubation time and temperature, and kit age can cause variation in binding.

12. This method can effectively eliminate the interference of the soluble receptors, binding proteins and other factors in biological samples.

# SAMPLE COLLECTION AND STORAGE

- Cell Culture Supernates Remove particulates by centrifugation.
- 2. **Serum** Use a serum separator tube (SST) and allow samples to clot for 30 minutes before centrifugation for 15 minutes at approximately 1000 x g. Remove serum, avoid hemolysis and high blood lipid samples.
- 3. **Plasma** Recommended EDTA as an anticoagulant in plasma. Centrifuge for 15 minutes at 1000 x g within 30 minutes of collection.
- 4. Assay immediately or aliquot and store samples at -20℃. Avoid repeated freeze-thaw cycles.
- Dilute samples at the appropriate multiple (recommended to do pre-test to determine the dilution factor).

Note: The normal rat serum or plasma samples are suggested to make a 1:2 dilution.

#### REAGENT PREPARATION

- 1. Bring all reagents to room temperature before use.
- 2. Wash Buffer Dilute 10mL of Wash Buffer Concentrate into deionized or distilled water to prepare 200mL of Wash Buffer. If crystals have formed in the concentrate Wash Buffer, warm to room temperature and mix gently until the crystals have completely dissolved.
- 3. Standard Reconstitute the Standard with 1.0mL of Standard /sample Diluent. This reconstitution produces a stock solution of 1000 pg /mL. Allow the standard to sit for a minimum of 15 minutes with gentle agitation prior to making dilutions. Pipette 500µL of Standard/sample Diluent into the 500 pg/mL tube and the remaining tubes. Use the stock solution to produce a 2-fold dilution series (below). Mix each tube thoroughly and change pipette tips between each transfer. The 1000 pg/mL standard serves as the high standard. The Standard/ sample Diluent serves as the zero standard (0 pg/mL).

If you do not run out of re-melting standard, store it at - 20℃. Diluted standard shall not be reused.

4. Working solution of Biotin-Conjugate anti-rat IL-4 antibody: Make a 1:100 dilution of the concentrated Biotin-Conjugate solution with the Biotin-Conjugate antibody Diluent in a clean plastic tube.

The working solution should be used within one day after dilution.

Working solution of Streptavidin-HRP: Make a 1:100 dilution of the concentrated Streptavidin-HRP solution with the Streptavidin-HRP Diluent in a clean plastic tube.

The working solution should be used within one day after dilution.

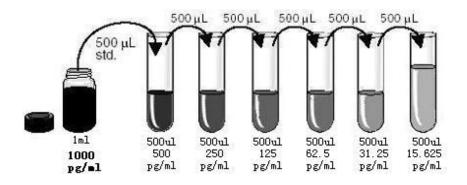


Figure 1: Preparation of IL-4 standard dilutions

#### **GENERAL ELISA PROTOCOL**

- 1. Prepare all reagents and working standards as directed in the previous sections.
- Determine the number of microwell strips required to test the desired number of samples plus appropriate number of wells needed for running blanks and standards. Remove extra microwell strips from holder and store in foil bag with the desiccant provided at 2-8°C sealed tightly.
- 3. Add 100μL of Standard, control, or sample, per well. Cover with the adhesive strip provided. Incubate for 1.5 hours at 37°C.
- 4. Aspirate each well and wash, repeating the process three times for a total of four washes. Wash by filling each well with Wash Buffer (350µL) using a squirt bottle, manifold dispenser or autowasher. Complete removal of liquid at each step is essential to good performance. After the last wash, remove any remaining Wash Buffer by aspirating or decanting. Invert the plate and blot it against clean paper towels.

- 5. Add 100 μL of the working solution of Biotin-Conjugate to each well. Cover with a new adhesive strip and incubate 1 hours at 37°C.
- 6. Repeat the aspiration/wash as in step 4.
- Add 100 μL of the working solution of Streptavidin-HRP to each well. Cover with a new adhesive strip and incubate for 30 minutes at 37°C Avoid placing the plate in direct light.
- 8. Repeat the aspiration/wash as in step 4.
- 9. Add 100 μL of Substrate Solution to each well. Incubate for 10-20 minutes at 37°C. Avoid placing the plate in direct light.
- 10. Add 100  $\mu$ L of Stop Solution to each well. Gently tap the plate to ensure thorough mixing.
- 11. Determine the optical density of each well immediately, using a microplate reader set to 450 nm.(optionally 630nm as the reference wave length;610-650nm is acceptable)

# ASSAY PROCEDURE SUMMARY

Prepare all reagents and standards as directed



Add 100µl standard or samples to each well, incubate 90 minutes,37℃

Add 100μl working solution of Biotin-Conjugate anti-rat IL-4 antibody to each well, incubate 60 minutes,37°C

∏ Aspirate and wash 4 times

Add 100µl working solution of Streptavidin-HRP to each well, incubate 30 minutes.37°C

Aspirate and wash 4 times

Add 100µl Substrate solution to each well, incubate 10-20 minutes,37℃.Protect from light.



Add 100µl Stop solution to each well. Read at 450nm within 30 minutes.

Figure 2: Assay procedure summary

#### **TECHNICAL HINTS**

- When mixing or reconstituting protein solutions, always avoid foaming.
- To avoid cross-contamination, change pipette tips between additions of each standard level, between sample additions, and between reagent additions. Also, use separate reservoirs for each reagent.
- 4. To ensure accurate results, proper adhesion of plate sealers during incubation steps is necessary. Substrate Solution should remain colorless until added to the plate. Stop Solution should be added to the plate in the same order as the Substrate Solution. Keep Substrate Solution protected from light. Substrate Solution should change from colorless to gradations of blue.
- A standard curve should be generated for each set of samples assayed. According to the content of tested factors in the sample, appropriate diluted or concentrated samples, it is best to do pre-experiment.

#### **CALCULATION OF RESULTS**

- Average the duplicate readings for each standard, control, and sample and subtract the average zero standard optical density.
- 2. Create a standard curve by reducing the data using computer software capable of generating a four parameter logistic (4-PL) curve-fit. As an alternative, construct a standard curve by plotting the mean absorbance for each standard on the y-axis against the concentration on the x-axis and draw a best fit curve through the points on the graph.
- 3. The data may be linearized by plotting the log of the IL-4 concentrations versus the log of the O.D. and the best fit line can be determined by regression analysis. This procedure will produce an adequate but less precise fit of the data. If samples have been diluted, the concentration read from the standard curve must be multiplied by the dilution factor.
- This standard curve is provided for demonstration only. A standard curve should be generated for each set of samples assayed.

Table 2: Typical data using the IL-4 ELISA (Measuring wavelength: 4 0nm, Reference wavelength: 630nm)

Standared(pg/ml)	OD.	OD.	Average	Corrected
0	0.024	0.022	0.023	
15.625	0.209	0.206	0.208	0.218
31.25	0.283	0.278	0.281	0.279
62.5	0.411	0.406	0.409	0.399
125	0.669	0.662	0.666	0.625
250	1.047	1.040	1.044	1.024
500	1.573	1.576	1.575	1.611
1000	1.946	1.949	1.948	1.941

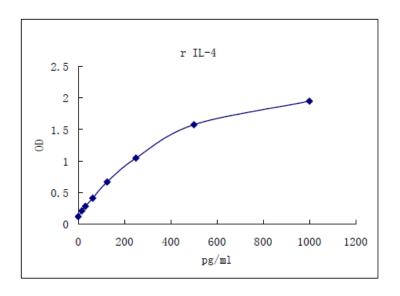


Figure 3: Representative standard curve for IL-4 ELISA. IL-4 was diluted in serial two-fold steps in Sample Diluent.

Do not use this standard curve to derive test results. A standard curve must be run for each group of microwell strips assayed.

#### BACKGROUND

Interleukin 4 (IL-4) is a pleiotropic cytokine produced primarily by activated T lymphocytes, mast cells and basophils (1-3).

The cDNA sequence of rat IL-4 predicts a 147 amino acid (aa) residue precursor protein containing a 24 aa residue signal peptide that is cleaved to generate the 123 aa residue mature protein containing four potential N-linked glycosylation sites and seven cysteine residues (4-6). Mature rat IL-4 shares 60% and 44% amino acid sequence similarity with mouse and human IL-4, respectively (5,6).

IL-4 has multiple immune response-modulating activities on a variety of cell types. It is an important regulator of isotype switching, inducing IgE production in B lymphocytes. It is an important modulator of the differentiation of precursor T helper cells to the Th2 subset that mediates humoral immunity and modulates antibody production. In addition, IL-4 has also been shown to have anti-tumor activity both in vivo and in vitro (1-3).

The biological effects of IL-4 are mediated by specific cell surface receptor complexes. Although IL-4 R does not bind IL-13 directly, it has been shown to complex with the low-affinity IL-13 R to form the functional high-affinity receptor complex for IL-13 (7, 8). In addition to the membrane-bound form of IL-4 R, a naturally occurring soluble form of IL-4 R has been identified in human and mouse biological fluids and in mouse cell culture supernates (9-11). Soluble IL-4 R has been to shown to bind IL-4 with high affinity in solution.

#### REFERENCES

- 1. Howard, M. and N. Harada (1994) in Guidebook to Cytokines and Their Receptors, Nicola N.A. ed., Oxford University Press, New York, p. 44.
- 2. Banchereau, J. and M.E. Rybak (1994) in The Cytokine Handbook, 2nd ed., A. Thomson ed., Academic Press, New York, p. 99.
- 3. Yokota, T. et al. (1990) in Peptide Growth Factors and Their Receptors I, Sporn, M.B. And A.B. Roberts eds., Springer-Verlag, New York, p. 577.
- 4. Carr, C. et.al., (1991) Biochemistry 30:1515.
- 5. McKnight, A.J. et.al., (1991) Eur. J. Immunol. 21:1187.
- 6. Richter, G. et.al., (1990) Cytokine 2:221.
- 7. Lin, J-X. et al. (1995) Immunity 2:331.
- 8. Zurawski, S.M. et al. (1995) J. Biol. Chem. 270:13869.
- 9. Mosley, B. et al. (1989) Cell 59:335.
- 10. Fernandez-Botran, R. and E.S. Vitetta (1990) Proc. Natl. Acad. Sci. USA 87:535.
- 11. Fanslow W.C. et al. (1990) Cytokine 2:3028.